


TDXPress PEBAX BALLOONS

Part Number	Material	Dia. (A) (mm)	Balloon Length (B) (mm)	Sleeve Length (mm)	Prox/Distal ID (inches) +/- .002"	Prox/Distal Cones (degrees)	Rated/Avg. Burst (atm)
40130-02	Pebax	2.0	8	5	.035/.035	20/20	12/16
40130-12	Pebax	2.0	20	5	.035/.035	20/20	12/16
40130-13	Pebax	2.0	30	5	.035/.035	20/20	12/16
40130-06	Pebax	3.0	15	5	.042/.042	20/20	12/16
40130-14	Pebax	3.0	20	5	.042/.042	20/20	12/16
40130-15	Pebax	3.0	40	5	.042/.042	20/20	12/16
40126-09	Pebax	4.0	20	5	.092/.092	34/34	10/14
40126-08	Pebax	4.0	30	5	.092/.092	34/34	10/14
40126-01	Pebax	4.0	40	5	.092/.092	34/34	10/14

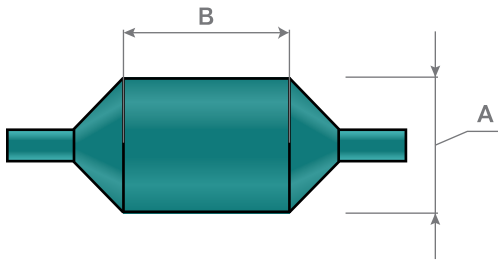

Pebax is a non-compliant balloon material that:

- can be used for medium to high pressures, but will require a thicker wall than PET
- tends to be more durable
- is a clear material
- is typically sterilized by means of Eto
- can be thermally welded or adhesive bonded

9/8/10

TechDevice Standard Products – Short Lead Times
CORPORATION



ISO 13485
CERTIFIED

1-888-TECHDEV • www.techdevice.com • 617-972-5800

FDA
REGISTERED