


# TDXPress URETHANE BALLOONS

Part Number	Material	Dia. (A) (mm)	Balloon Length (B) (mm)	Sleeve Length (mm)	Prox/Distal ID (inches) +/- .002"	Prox/Distal Cones (degrees)	Rated/Avg. Burst (atm)
40131-02	Urethane	2.0	8	5	.035/.035	20/20	12/16
40131-12	Urethane	2.0	20	5	.035/.035	20/20	12/16
40131-13	Urethane	2.0	30	5	.035/.035	20/20	12/16
40131-06	Urethane	3.0	15	5	.044/.044	20/20	12/16
40131-14	Urethane	3.0	20	5	.044/.044	20/20	12/16
40131-15	Urethane	3.0	40	5	.044/.044	20/20	12/16
40127-09	Urethane	4.0	20	5	.093/.093	34/34	10/14
40127-08	Urethane	4.0	30	5	.093/.093	34/34	10/14
40127-01	Urethane	4.0	40	5	.093/.093	34/34	10/14

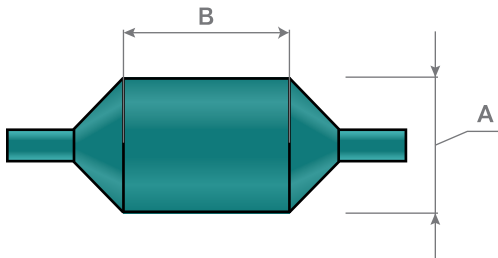

Urethane is a semi-compliant balloon material that:

- is used for lower pressure applications
- is typically sterilized by means of Eto
- can be thermally welded or adhesive bonded

3/25/11

**TechDevice** Standard Products – Short Lead Times  
CORPORATION



ISO 13485  
CERTIFIED

1-888-TECHDEV • www.techdevice.com • 617-972-5800

FDA  
REGISTERED